

TrioMed Antimicrobial Medical Adhesive Tape

INTENDED PURPOSE

The TrioMed Antimicrobial Medical Adhesive Tape is a single use disposable device intended for general securement of wound dressings, tubing or any other device. The tape is treated with a patented iodine-based antimicrobial to protect the external surface from microbial deterioration. Hypoallergenic.

INSTRUCTIONS:

Product can be worn for up to 5 days, or as per instructions from your nurse/physician.

Product has a 5 year shelf life.

For external use only.

Single use. Once removed it must be discarded as it will not adhere appropriately to skin anymore.

Apply on clean, dry and oil or cream free skin.

Clean and dry the skin over intended application area.

Unwind the desired length of medical tape to maintain adhesion, avoid touching the adhesive area of the tape as much as possible.

Customize the tape size and shape as required for your application.

Start by applying one side of the tape to the patient's skin. Without stretching or pulling, smooth down the remaining segment of tape and ensure it has adequately adhered with no wrinkles or gaps in adhesion.

Store at room temperature away from direct sunlight.

Dispose of according to your country's regulations.

CONTRAINDICATIONS:

Discontinue use and consult a physician if irritation or sensitivity develops.

This product contains iodine, it should not be used by persons with known iodine hypersensitivity.

WARNINGS:


The product helps protect against certain particulate contaminants but does not eliminate exposure to or risk of contracting the disease or infection.

Any serious incident that occurs when using this device should be reported to the manufacturer and the national competent authority.

PRODUCT BENEFITS

The outer active layer with TrioMed antimicrobial technology (iodine based) is to protect the external surface from microbial deterioration. The iodine is embedded in the external surface therefore does not leach or escape the material, its purpose is to render bacteria and viruses inactive when they come into contact.

Explanation of symbols:

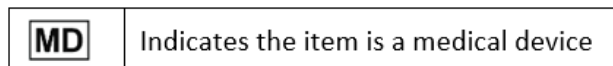

	 	Obelis S.A. Bd. Général Wahis 53 1030 Brussels, Belgium Tel: +(32) 2 732-59-54 Fax: +(32) 2 732-60-03 mail@obelis.net	 	I3 BioMedical Inc. 14163 boul. Curé Labelle Suite 50 Mirabel, Québec, Canada J7J 1M3 WWW.I3BIOMEDICAL.COM
---	--	--	--	---